




Life, we care
Health, we serve
We reach everywhere



Guangdong Shantou Minston Medical Instruments Co., Ltd.

Electric operating table MST-JH99

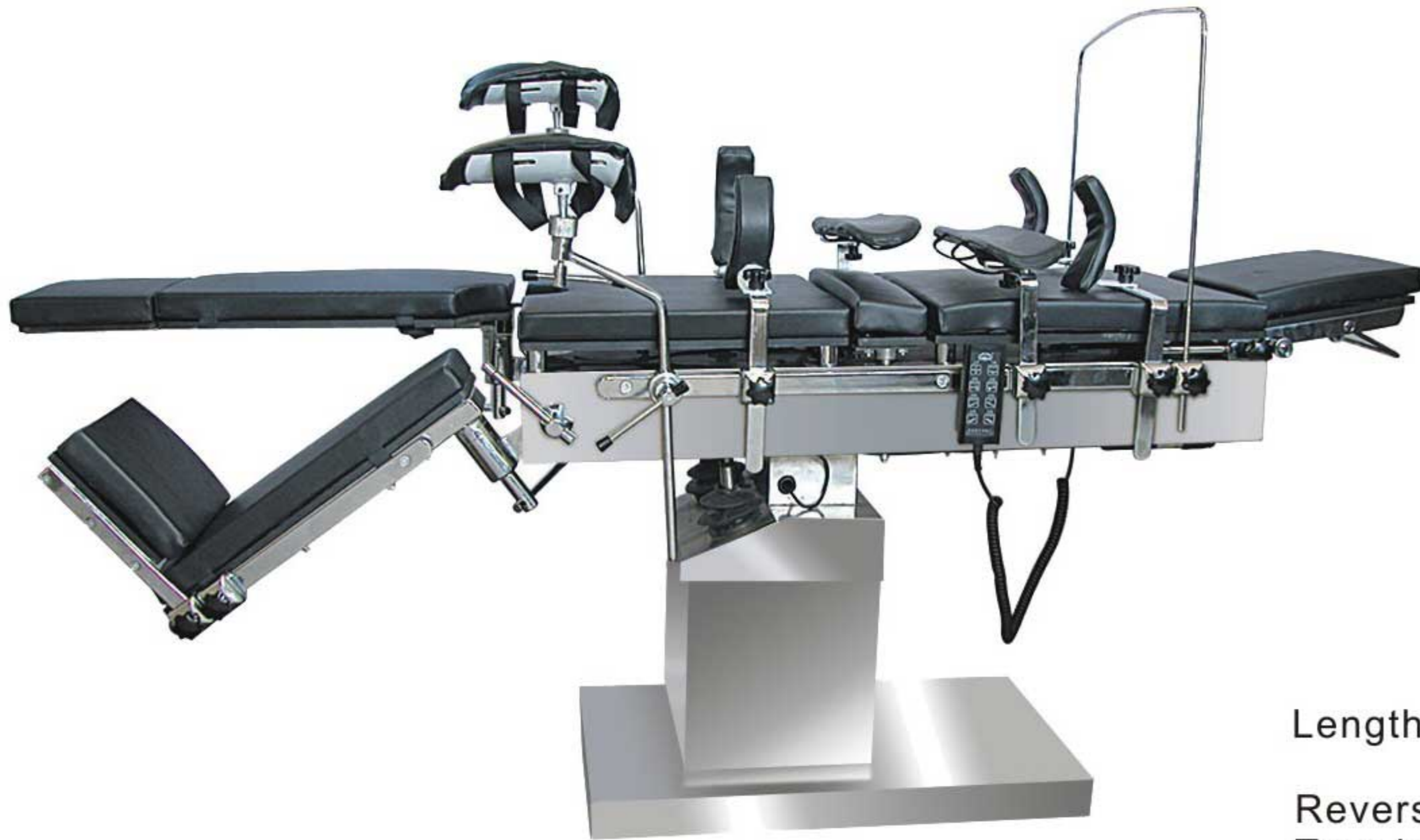

Electric Operating Table is designed for all kinds of operations in hospitals. Besides all functions of the Chinese-built operating tables. It has its own advantages:

1. Leg section can be rotated, opened, disassembled and adjusted easily. It is most suitable for the urological operation.
2. The table top, made of high strength plastic plate, can move backwards and forwards, X-ray examination can be carried out on it.
3. With 24V direct supply, the table is safe and reliable. The table's elevation, lowering, reversed trendelenburg; trendelenburg, lateral tilt are all automatically controlled by the Button.
4. This table can be disposed different kinds of table tops for different use. Electric Operating Table is only such kind of product in China. It is one of the most advanced operating table with excellent functions.

Length:2100mm Width:480mm Max Height:1000mm Min Height:750mm


Reversed Trendelenburg:≥25°
Trendelenburg:≥25°
Lateral Tilt:≥15°
Leg section Opened ≥90°
Head Section lowered from the Horizontal ≥90°

Leg section Lowered from the Horizontal:≥90°
Kidney Bridge Elevation:≥60mm
Backward and Forward Movement:>400mm
Head Section Raised from the Horizontal ≥45°
Supply AC 220V±10%,50HZ

Electric operating Table MST-JH2000G

Electric operating Table was exploited independently by our company with innovative and advanced technologies introduced from domestic and overseas. This table can be used in almost all kinds of operation including orthopedics, otology, gynaecology, urology etc. It is a world leading product.

This table is composed by electric Hydraulic pressure system, Control circuit and X-ray transmissive head section. Full function, easy adjustment, with the help of C-arm, the table can be more capable for any operation. The back section is made form five parts, which can be easily adjusted for different operations such as neck, lumbar and spine operations. With the good performance, you will find it is a perfect operation equipment.

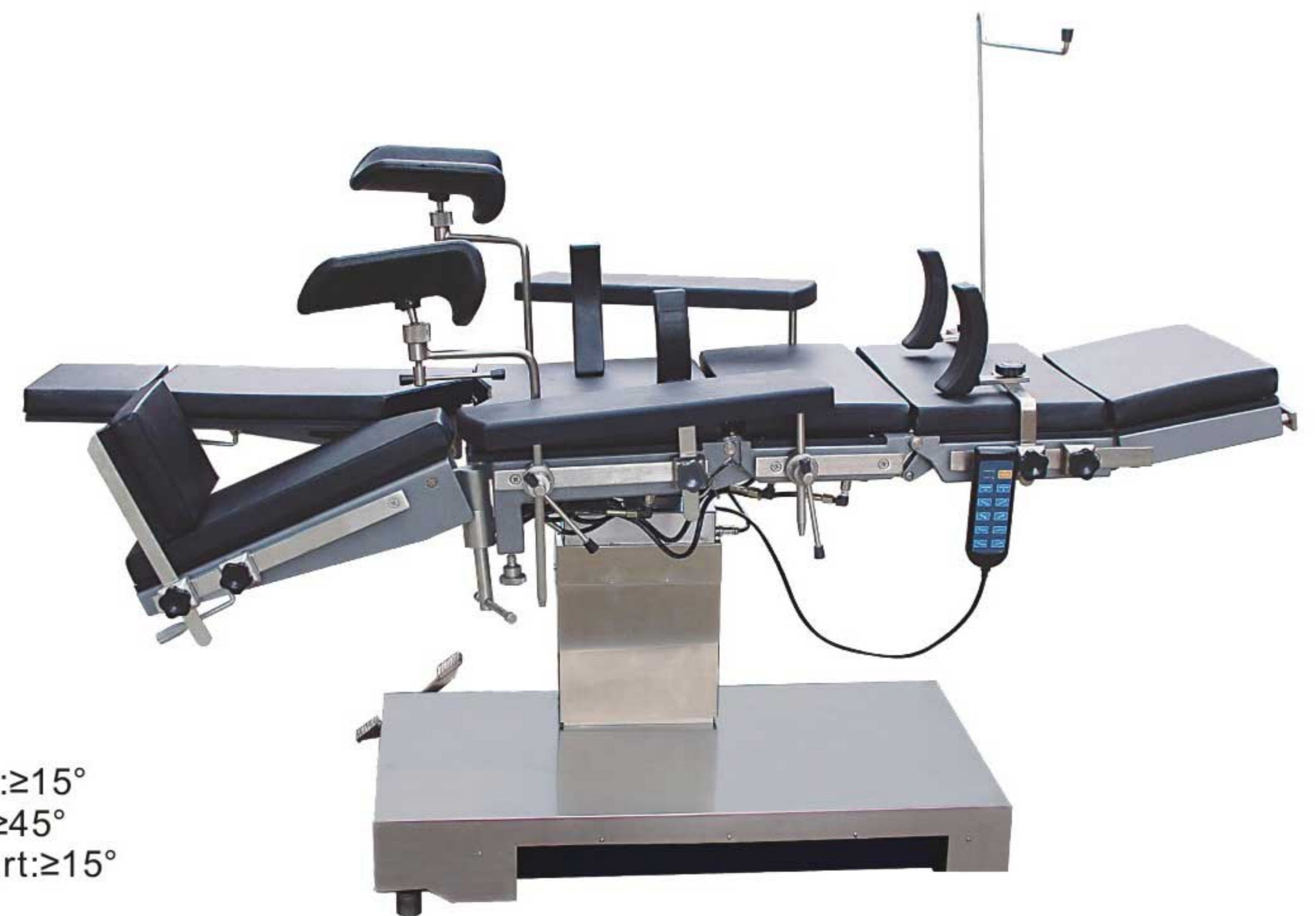


Length:2000mm Width:500mm Max Height:750~1000mm Adjust:500mm

Reversed Trendelenburg:≥25°
Trendelenburg :≥25°
Lateral Tilt:≥20°
Folds Upward of Head Section:≥35°
Folds Downward of Head Section:≥90°
Folds Upward of back Section's Foreside:≥45°

Folds Downward of back Section's Foreside:≥15°
Folds Upward of back Section's Lower part:≥45°
Folds Downward of back Section's Lower part:≥15°
Folds Downward of Leg section:≥90°
Eversion Leg section:≥90°
Supply AC 220V±10%,50HZ

Multi-purpose operating table, Side-controlled MST-JH3008AB

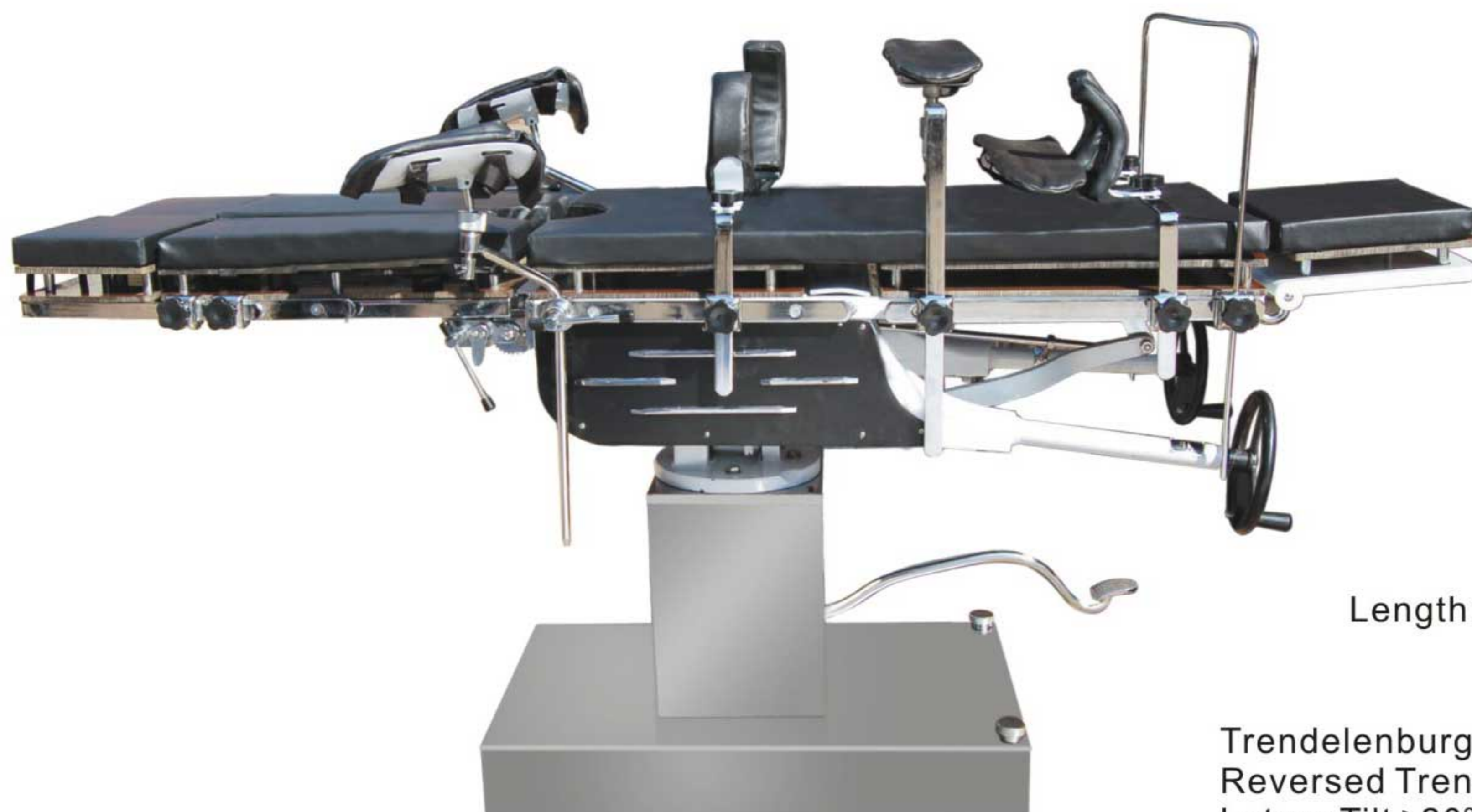

This table is designed for surgical procedures of the head, neck, thorax and abdomen, and extremities as well as for gynaecological, oto-thinology and orthopaedics operations. It can be raised or lowered hydraulically by means of an oil pump. Various positions can be controlled on the head. Model 3008AB made of double decker plastic, X-ray examination can be available during the operation. Leg section can be opened 90 and disassemble. It is most for the urological operation, It was awarded "High Quality Product" certificated by Jiangsu Province. Chassis inside & outside cover of model 3008AB are stainless steel.

Length:2040mm Width:480mm Height:800~1040mm

Trendelenburg ≥30°
Reversed Trendelenburg ≥50°
Lateral Tilt ≥20°
Back Section Raised From the Horizontal ≥80°
Back Section lowered From the Horizontal ≥8°

Head Section Raised from the Horizontal ≥30°
Head Section lowered from the Horizontal ≥30°
Leg section Lowered from the Horizontal ≥90°
Kidney Bridge Elevation ≥120mm
Leg section Opened ≥50°
Lowered from the Horizontal ≥90°

Multi-purpose operating table, Side-controlled MST-JH3008B



This table is designed for surgical procedures of the head, neck, thorax and abdomen, and extremities as well as for gynaecological, oto-rhinolaryngological and orthopaedics operations. It can be raised or lowered hydraulically by means of an oil pump. Various positions can be controlled on the head. Leg section can be opened 90° and disassemble. It is most for the urological operation, It was awarded "High Quality Product" certificated by Jiangsu Province. Chassis inside & outside cover of model 3008B are stainless steel.

Trendelenburg $\geq 30^\circ$
 Reversed Trendelenburg $\geq 50^\circ$
 Latera Tilt $\geq 20^\circ$
 Back Section Raised From the Horizontal $\geq 80^\circ$
 Back Section lowered From the Horizontal $\geq 8^\circ$

Head Section Raised from the Horizontal $\geq 30^\circ$
 Head Section lowered from the Horizontal $\geq 90^\circ$
 Leg section Lowered from the Horizontal $\geq 90^\circ$
 Kidney Bridge Elevation $\geq 120\text{mm}$

Length: 2040mm Height: 800~1040mm Width: 480mm

Multi-purpose operating table, Side-controlled MST-JH3008D

Character:
 Suit for all surgery operating or examing
 It take hydraulic pump by footboard
 It take cylinder to control each action, easy to operate and localize
 New configuration is more sound
 Bulit in wasit section is easy to use
 Mattress is manufactured with anti-ccid and alkali material
 3008D is one multi function operating table for thoracal, ophthalmology, urinology operation ect, It is a good all-round equipment for clinical use.



Head Section Raised from the Horizontal $0^\circ \sim 35^\circ$ (Detachable)
 Head Section lowered from the Horizontal $0^\circ \sim 40^\circ$ (Detachable)
 Back Section Raised From the Horizontal $0^\circ \sim 70^\circ$
 Back Section Raised From the Horizontal $0^\circ \sim -10^\circ$
 Leg section Lowered from the Horizontal $0^\circ \sim 75^\circ$ (Detachable)
 Reversed Trendelenburg $0^\circ - 15^\circ$
 Trendelenburg $0^\circ \sim 15^\circ$
 Left $0^\circ - 10^\circ$
 right $0^\circ - 10^\circ$
 Adjust of Waist section $\geq 80\text{mm}$
 Load weight 135Kg
 Baseplate removed or locked

Length: 1990mm Width: 480mm Min Max Height: $\geq 940\text{mm}$ Height: $\leq 700\text{mm}$ Adjust: $\geq 240\text{mm}$

Ordinary operating table MST-JHBS



This table is suitable for surgical procedures of the head, neck, thorax and abdomen, and extremities as well as for gynaecological, oto-rhino-laryn-gological and orthopaedics operations. It is characterized by versatility, easy operation, good and cheap. Chassis and inside & outside cover of model BS are stainless steel.

Length: 2000mm Height: 750~1000mm Width: 480mm

Back Section Raised From the Horizontal $\geq 45^\circ$
 Back Section lowered From the Horizontal $\geq 20^\circ$
 Head Section Raised from the Horizontal $\geq 45^\circ$

Head Section lowered from the Horizontal $\geq 90^\circ$
 Leg section Lowered from the Horizontal $\geq 90^\circ$
 Kidney Bridge Elevation $\geq 100\text{mm}$